

T2 mk2 Operation - a first look

Target Meeting 24th Feb 2010
Paul Hodgson

BCD vs Time Total

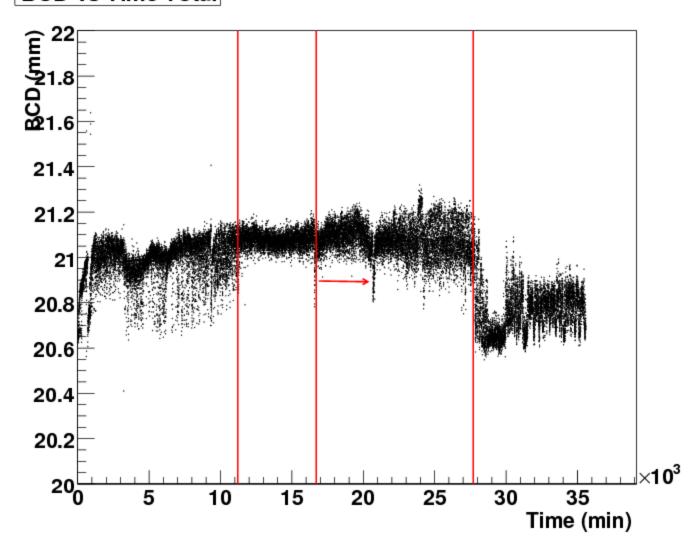

Acceleration vs Time Total

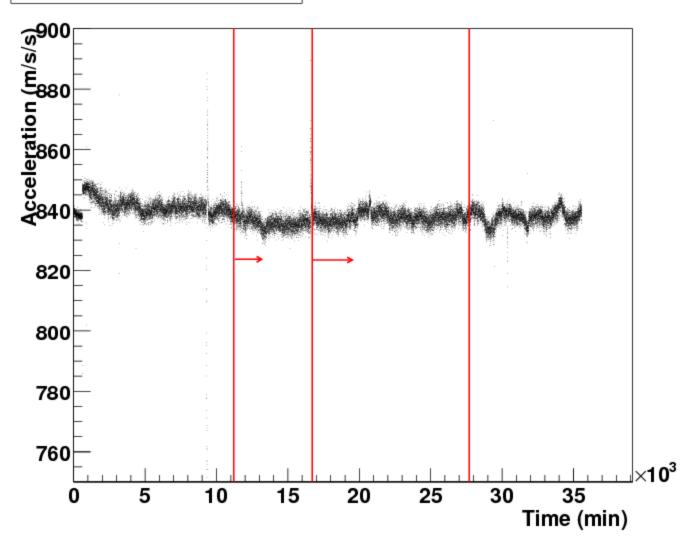

- See 4 possibly 5 distinct regions
- 11200 min 670 k pulses
- Target bedding in ?
- 16700 min 1×10⁶ pulses
- 20000 min 1.2×10⁶ pulses
- Good stable operation
- 27700 min 1.66×10⁶ pulses
- Increasing BDC rms (width) until rapid change
- Material being displaced after prolonged wear
- To end 2.15×10⁶ pulses
- Decrease in mean BCD, was 21.1mm, reduced to 20.8mm

